

Cue[®] Cell Processing System

Unmatched precision and control make Cue the definitive instrument for small volume cell processing.



Final product container(s) are user provided and are not part of the Cue consumables.
Cue Light Shield not shown here.

For Laboratory Use only

Scan Here or
Visit [choosecue.com](https://www.choosecue.com)
to Learn More



Cue reflects our commitment to developing tools that meet your research needs.



Cryopreservation Preparation¹

Concentrate and wash cells into a buffer. Then dilute with cryoprotectant, achieving desired cell concentration, and generate dose and sample aliquots in freezing bags.



Small Volume Culture Media Exchange

Concentrate and wash cells into fresh media for additional culture and dispense aliquots directly into multiple culture bags or vessels.



Electroporation Preparation


Concentrate and wash cells into electroporation buffer to achieve desired cell concentration for electroporation.

Novel technology to solve the most challenging problems

Cue's functionally closed syringe pumps set a new standard of precision and control for small volume cell processing. By pairing this breakthrough with Lovo's proven spinning membrane filtration technology, we can now tackle some of the most important steps in a manufacturing process. Whether you're formulating and aliquoting a final product or performing upstream wash and concentration, Cue delivers unmatched results across a wide range of critical applications.



Instrument Components

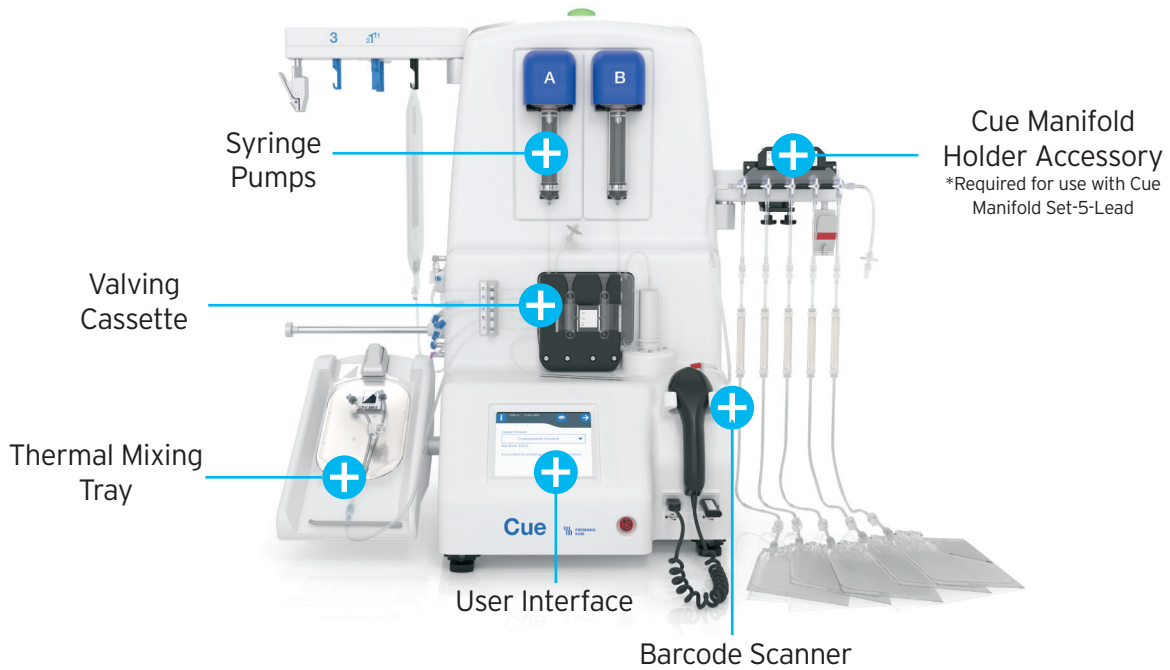

Take control of processing consistency and product quality

The combination of our functionally-closed syringe pumps, a valving cassette, and our spinning membrane separation technology enable high accuracy concentration, formulation, and aliquoting. Additionally, Cue's highly configurable protocol design suite gives users the ability to customize protocols and determine final and in-process volumes required before procedure start to achieve concentration and formulation ratio targets.

Building blocks make protocol design easier than ever

Cue performs a wide range of cell processing functions. With the included Cue Desktop Application Software, you can easily create protocols that meet your processing needs. You can choose to do all steps, or only a few. Please speak with your local Application Specialist to learn more.

- Concentrate
- Wash
- Dilute
- Formulate
- Aliquot

Viable Cell Recovery ¹ 98.3 ± 4.9	Automated Source Processing Time ¹ 22.0 ± 1.8 min
Viability Change (from Source) ¹ -0.6 ± 2.1%	Aliquot Concentration Variability ¹ 10.8 ± 4.6%

Cue Processing System

Instrument features

Volumes

Source Input	Validated: 10 - 500 mL (Software entry allows 1 - 9,999 mL)*
Minimum Final Product Volume	10mL (Single Aliquot), 2 mL (Multiple Aliquots)
Volume Accuracy	±10% or 1mL, whichever is greater

*The system allows for entry of up to 9,999 mL however ability to process may depend on cell type, cell concentration, and source volume

General Features

- Active Temperature Control Range: 3-22 °C
- Protocol Options*: Concentrate, Wash, Dilute, Formulate, Aliquot
- Save up to 20 protocols on the instrument at any time

*Protocols can be configured to include all or a subset of the available options

Cue Desktop Application Software Included With Purchase

- Protocol design and management
- Offline procedure estimator
- View and store procedure records

System Components

94 lb, 35.3" X 18.8" X 32.6" (W x D x H*) benchtop instrument
42.6 kg, 89.7 X 47.8 X 82.8 cm (W x D x H*) benchtop instrument

*The instrument height may be variable depending on the configuration of the instrument feet. The maximum height is listed above.

References

¹Data on File at Fresenius Kabi USA

Product Ordering Information

Product	Product Description	Quantity/Case
X6R5004	Cue Primary Set - 4µm	4
X6R5002	Cue Manifold Set - 8-Lead	10
*X6R5005	Cue Manifold Set 5-Lead	4
X6R5006	Cue Sampling Set	6
*6112241039S	Cue Manifold Holder	1

*Cue Manifold Set 5-Lead Requires the Cue Manifold Holder that is purchased separately in order to use this product.

Not all products listed or shown are available or approved for sale in all countries. Contact your local representative for additional information and product availability.

The Cue Cell Processing System is for laboratory use only and may not be used for direct transfusion. Appropriate regulatory clearance is required by the user for clinical use.

Refer to the Cue Cell Processing System User's Guide for a complete list of warnings and precautions associated with the use of these products.

**Learn how choosing Cue can help make your lab more productive:
Call 1.800.333.6925 or [choosecue.com](https://www.choosecue.com)**

Trademarks referred to are property of their respective owners.



**FRESENIUS
KABI**

US Headquarters:
Fresenius Kabi
Three Corporate Drive
Lake Zurich, IL 60047