



**FRESENIUS  
KABI**

caring for life



## CompoGuard®

### Whole blood donation mixing/weighing device

Donation, quality and safety combined in whole blood processing that includes increased data management flexibility with digital enhanced cordless technology (DECT).

High capacity Lithium battery technology for increased donation capabilities

Flexible, mobile design with comfortable workstation

Streamlined data management with DECT supported bi-directional data transfer

# Trusted features of our

## Optenna (optical alarm extension) for optimal donation supervision

- The donation event can be monitored from a distance
- The blinking optenna leads you directly to the correct donor when action is needed

## Donor engagement to improve donation process


- Donors can directly monitor their donation progress

## In-line tray for optimal results

- Integrated filter holders to support proper mixing
- Large tray for optimal bag positioning and accurate weighing

## Safety clamp prevents air from entering the collection bag

- The donation line tube is fixed inside the safety clamp
- The tube clamp can only be opened after successful donation and sealing
- The heat sealing is integrated and documented by the CompoGuard

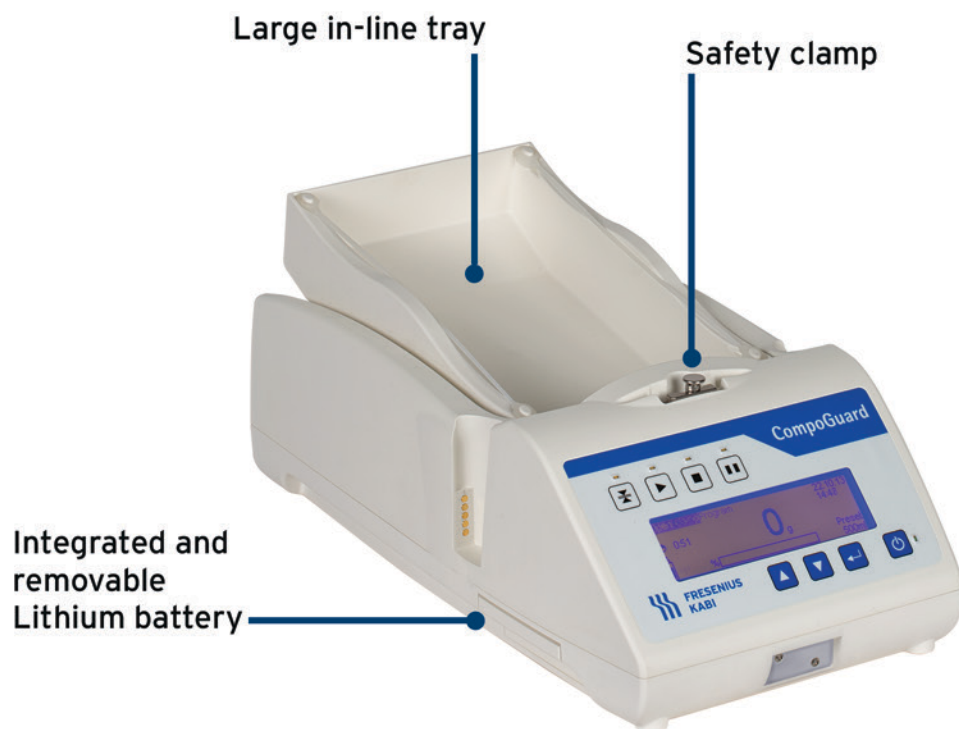

## CompoGuard® Basic (100 - 240 V)

- Large in-line tray
- Safety clamp
- Integrated Lithium battery

# CompoGuard models



## CompoGuard® Data (100 - 240 V)

Basic CompoGuard device plus:

- Reliable and safe bi-directional data transfer
- Ergonomic working station
- DECT option
- Gooseneck
- USB port



## CompoGuard® Complete (100 - 240 V)

Our top model combines CompoGuard Data plus:

- Reliable and safe bi-directional data transfer
- Ergonomic working station
- Reliable sealing function and documentation
- Rechargeable battery

### Flexible data handling with DECT technology

- Bi-directional data transfer
- Easy and flexible programming
- Reliable data recording, surveillance and documentation

In addition to DECT and cable transmission, data transfer via USB stick is possible



# CompoGuard

More donations.

More flexibility.

More benefits...

## Increased number of donations per fully charged battery

- New high capacity Lithium battery technology
- Enhanced capability of approximately 70 donations per battery charge for mobile donations
- External charging station available
- Battery easily exchangeable



## Increased data management flexibility

- Donation data now able to be downloaded via dedicated Fresenius Kabi USB stick



## One further step towards paperless donation

- A new software feature enables the direct donor verification via name, surname and date of birth

## Mobile donation system

- Robust transport case for mobile donation

## Comfortable working station

- The Gooseneck allows the operation of the CompoGuard without bending down
- The large in-line tray and the tube clamp enable comfortable loading of the blood bag
- Transport case is also mobile and comfortable table during donation



# CompoGuard® Whole blood donation

---

## Ordering Information

For more information such as literature, technical details and working procedures, please contact your local sales representative.

Process efficiency to  
help you achieve more

CE marked, Class IIa according Directive 93/42/EEC, as amended 2007/47/EC.

Electrical safety, IEC 60601-1

Cleared for sales in the United States

The signs/names marked with ® are registered trademarks of the Fresenius Group in selected countries.



**FRESENIUS  
KABI**

caring for life

Fresenius Kabi AG  
61346 Bad Homburg  
Germany

Phone: +49 6172 608-0  
Fax: +49 6172 608-5786

compoguard@fresenius-kabi.com  
www.fresenius-kabi.com