MEDICAL DEVICES

CompoLab TM

FAST, ERGONOMIC AND RELIABLE HEMOGLOBIN ANALYZER







CompoLab TM

Small handheld hemoglobin analyzer with less than 2 seconds measuring time

Fast, Ergonomic and Easy to Use

- Display always in stand-by mode device ready to use
- No moving parts robust for mobile use
- Light weight approx. 180 g
- No calibration needed minimal maintenance and cleaning
- Optional Fresenius Kabi software "HB-Master" can be installed on local PC - easy data transfer via USB cable
- Optional barcode scans simply connect a barcode reader to the PC
- 2 charging possibilities USB and adapter plug




Fast Measurement within 2 sec.

Small Size with Large Display

Long Lasting Battery Capacity

Optional Data Management



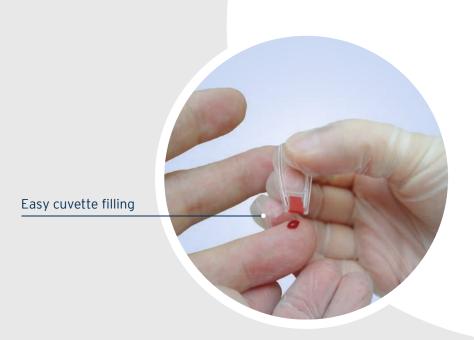

Reagent free cuvette with unique design and long expiration dating

- No chemical reaction inside the cuvette -Insensitive to temperature and moisture
- Air vent prevents formation of bubbles in the cuvette
- Cuvette shelf life 2.5 years even in opened foil bag
- Can be filled from any angle
- Only one drop of blood needed
- Same cuvettes for CompoLab TM and CompoLab TS







Ordering information

For more information such as literature, technical details, working instructions, as well as equipment, please contact your local sales representative or call us directly.

ORDER NUMBERS

9034000 CompoLab TM

9033310 CompoLab TM/TS cuvettes, 100 cuvettes in foil bag (5 bags per case)

9034010 CompoLab TM transport case

TECHNICAL DATA

Measuring range 0 to 25.5 g/dL (0 to 255 g/L; 0 to 15.9 mmol/L)

Measurement accuracy 1.8 - 25.4 g/dL within +/- 1.3% to ICSH

Storage temperature CompoLab TM: 0°C to +50°C (-30°C to +70°C during transport, 24h max.)

Cuvettes: 0°C to +50°C (-30°C to +70°C during transport, 24h max.)

Duration of use Several months, with a fully charged battery

Power supply 100-240 V AC, 50-60 Hz Dimensions 9 x 15 x 4 cm (W x D x H)

Weight Approx. 180 g

EC-Directive 98/79/EC

CE marked, IEC 61010 compliant



Fresenius Kabi AG

Else-Kroener-Strasse 1 61352 Bad Homburg, Germany Phone: +49 (0) 61 72 608 4993 Fax: +49 (0) 61 72 686 4728 compolab@fresenius-kabi.com www.fresenius-kabi.com Your local sales organization: