


MEDICAL DEVICES

CompoSeal Mobilea II

THE PORTABLE TUBE SEALER FOR HIGH QUALITY SEALS

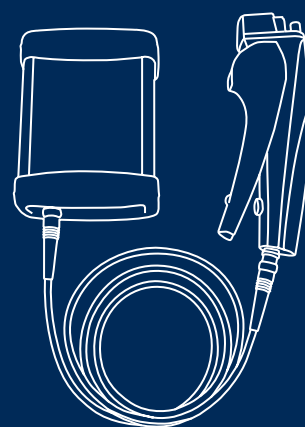

**FRESENIUS
KABI**

caring for life

CompoSeal Mobilea II

New technology for portable and safe sealing

Mobilea II meets the high demands of professional sealing and combines it with the benefits of modern battery technology. The Lithium Polymer battery realises a high sealing performance and offers a portable light weight solution. Mobilea II is therefore designed for mobile sealing in donation teams or at processing departments. Furthermore the intelligent power management combines high quality seals with short sealing times.

CompoSeal Mobilea II - the portable light weight sealer with a high battery capacity

- Modern Lithium Polymer battery provides a capacity of up to 1000 seals per full charge on standard blood bag tubes
- Designed for portable use: Low weight power pack of just 760g, robust metal housing, rubber pads for shock absorption
- High quality seals independent of the battery capacity
- Automatic adjustment of sealing time depending on tube thickness
- Approved for bedside sealing while donor is connected

CompoSeal Mobilea II Transport Case

The CompoSeal Mobilea is delivered in a durable and handy transport case.



Description	Dimensions (l x w x h)	Weight
CompoSeal Mobilea II set ¹⁾	395 mm x 106 mm x 295 mm	2364 g
CompoSeal Mobilea II powerpack	110 mm x 52 mm x 145 mm	760 g
CompoSeal Mobilea II handsealer	215 mm x 65 mm x 30 mm	315 g
CompoSeal Mobilea II power supply unit	65 mm x 75 mm ²⁾ x 130 mm	270 g
CompoSeal Mobilea II handsealer cable	1,470 mm x 5 mm x 5 mm ³⁾	83 g
Transport case for CompoSeal Mobilea II	395 mm x 106 mm x 295 mm	900 g

¹⁾ In transport case, ²⁾ incl. adapter, ³⁾ Length from plug to plug, Ø of the cable itself

- high quality seals
- reliable operation
- quick and easy handling
- process control ensures leakage detection



ORDERING INFORMATION

For more information such as literature, technical details, working instructions, as well as equipment, please contact your local sales representative or call us directly.

ORDER NUMBER CompoSeal Mobilea II

9027011 CompoSeal Mobilea II set

ORDER NUMBER Spare Parts

9027001 CompoSeal Mobilea II powerpack
9027021 CompoSeal Mobilea II handsealer
9027031 CompoSeal Mobilea II charger
9027041 CompoSeal Mobilea II hand sealer cable
9027051 CompoSeal Mobilea II transport case

Class I accord. EC-Directive 93/42/EEC
CE marked, IEC 60601-1 3rd Edition Compliant
Cleared for sales in the US



**FRESENIUS
KABI**

caring for life

Fresenius Kabi AG
Else-Kroener-Strasse 1
61352 Bad Homburg, Germany
Phone: +49 (0) 61 72 686 5047
Fax: +49 (0) 61 72 686 4728
composeal@fresenius-kabi.com
www.fresenius-kabi.com

Your local sales organization: