



CompoSeal® Slim Tube Sealing System

Flexible and efficient tube sealing technology

Versatile tube sealing platform for blood center, laboratory, and hospital use.

Cascade up to 8 devices for
segment sealing

Quick and easy to set-up, use and
re-configure

High-quality seals for reliable operations

Optional hand sealer adds flexibility
and efficiencies

CompoSeal Slim

Tube Sealing System

Flexible desktop tube sealing system for single and cascade sealing


Combine up to 8 CompoSeal Slim devices for cascade sealing

- Adjustable segment length (60 - 120 mm)
- Color coded LEDs guide the operator through the sealing process. The sealing process is designed to seal from the last segment toward the blood bag to prevent pressure build-up in the tube
- One power cord needed per 2 devices in cascade configuration



CompoSeal Slim - versatile platform for high quality tube sealing

- Seals tubing with ODs up to 5.9 mm
- Integrated tube positioning system with splash guard
- LEDs indicate sealing process
- Sealing time is approximately 1-2 seconds, depending on tube size
- Optional hand sealer for added flexibility and efficiency
- Sealing head cap detaches for easy cleaning and disinfection

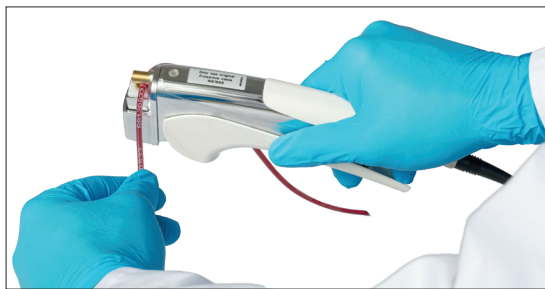

Range of sealable tubing 3.0-5.2 mm (outer diameter)

Optional hand sealer

- The ergonomic design of the hand sealer supports ambidextrous use
- Both low grip angle and low grip force activate the sealing process
- Visual and acoustic alerts guide the operator through every seal
- Optional splash guard offered for hand sealer

For more information, such as technical details, instructions for use, and training materials, please contact your local sales representative.

Refer to the Operating Instructions for a complete list of warnings and precautions associated with the use of this product.



**FRESENIUS
KABI**

caring for life

US Headquarters
Fresenius Kabi
Three Corporate Drive
Lake Zurich, IL 60047
800-333-6925

Trademarks referred to are property of their respective owners.

Copyright © 2023 Fresenius Kabi AG. All rights reserved.

0001-CSS-01-09/18
January 2023